


MEGADYNE™ Telescoping Smoke Evacuation Pencil

Optimized Device Performance Guide

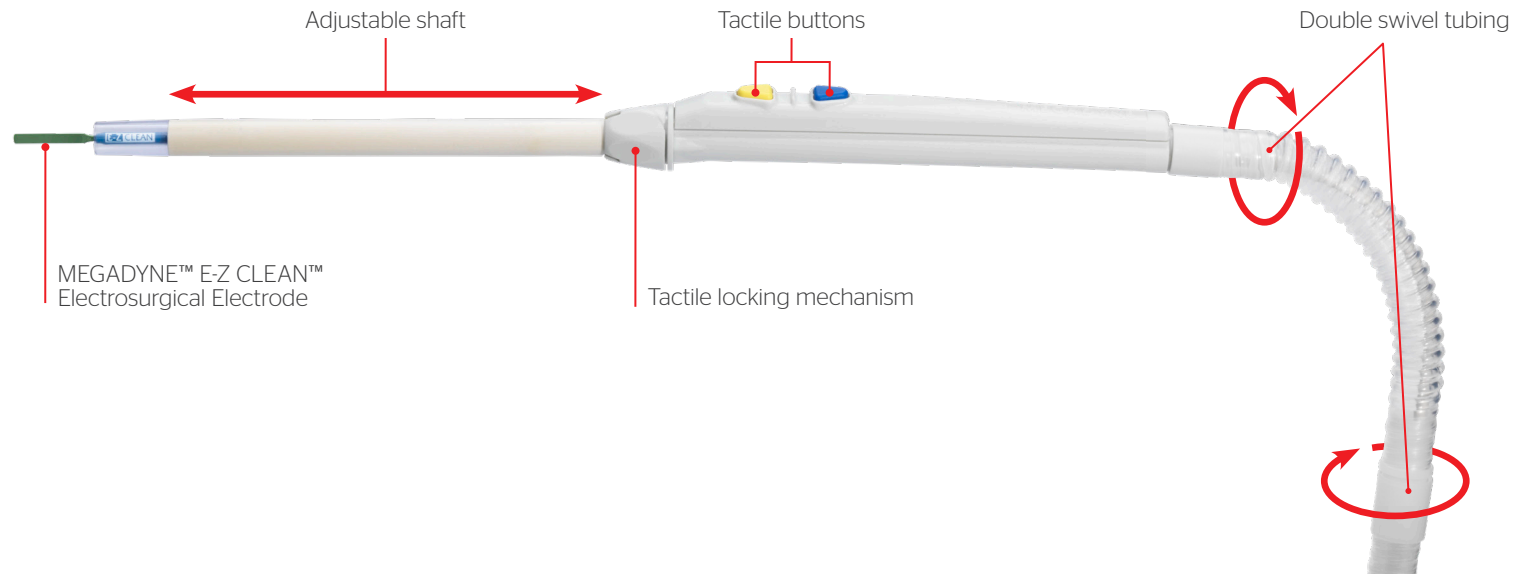

Before you begin

1. Check electrode and remove tip protector:

Ensure the electrode is fully inserted. Remove the tip protector from electrode and discard.



2. Attach to generator:

Insert the pencil plug into the receptacle on the MEGADYNE™ ElectroSurgical Generator.




3. Attach smoke tube to filter:

Securely connect the smoke evacuation tubing into the filter on the MEGADYNE™ Smoke Evacuator.



Shaft extension

- 1. Adjust the shaft:** Pull shaft away from the pencil shroud to extend to desired length.
- 2. Lock the locking mechanism:** After reaching desired shaft length, turn the locking mechanism to the left to lock.

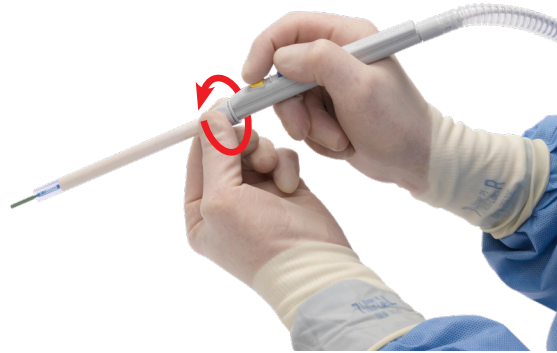

Clean the device

- 1.** MEGADYNE™ Telescoping Smoke Evacuation Pencil is coated with PTFE to reduce eschar buildup and to aid in the easy removal of eschar with a damp gauze or sponge.

CAUTION: Do not clean instrument with abrasives. Do not use a scratch pad or other abrasive cleaner to remove eschar. This could damage the PTFE coating.

