

Compression injury of the circular stapler for gastrointestinal end-to-end anastomosis: preliminary *in-vitro* study

Gyung Mo Son^{1,2}, Myeong Sook Kwon^{2,3}, Hong-min Ahn¹, In Young Lee², Gun Ho Kim⁴, Kyoung Won Nam⁴, JoonWon Lee⁵, Chang-Suk Kong³

¹Department of Surgery, Pusan National University School of Medicine, Yangsan, Korea

²Research Institute for Convergence of Biomedical Science and Technology, Pusan National University Yangsan Hospital, Yangsan, Korea

³Department of Food and Nutrition, College of Medical and Life Science, Silla University, Busan, Korea

⁴Department of Biomedical Engineering, Pusan National University School of Medicine, Yangsan, Korea

⁵Biomedical Materials Science, School of Dentistry, University of Birmingham, Birmingham, United Kingdom

Purpose: This preliminary *in-vitro* study was designed to evaluate the risk factors of compression injury from use of a circular stapler for end-to-end anastomosis.

Methods: Transparent collagen plates were prepared in dry and wet conditions. Physical properties of collagen plates and porcine colon tissue were examined using a rheometer. Adjustable and fixed-type circular staplers were applied on the collagen plates and the gap distance and compressive pressure were measured during anvil approximation. Tissue injury was evaluated using a compression injury scale. Compression properties were accessed to optimal or overcompression based on gap distance.

Results: Unacceptable injuries were rarely observed on the dry collagens, regardless of compression device. In the adjustable compression, the compressibility ratio was similar between dry and wet collagen. Overcompression and unacceptable injury increased on the wet collagens. In the fixed compression, the compressibility ratio increased significantly and unacceptable injuries were observed in more than 50% of wet collagens. Peak pressure was significantly higher in the fixed-compression types than those of adjustable type. On bivariate correlation analysis, fixed-compression type and wet collagens were respectively associated with overcompression. On multivariate analysis, edematous collagen condition was the most important risk factor and proximal anvil side, fixed compression type, and overcompression were also independent risk factors for unacceptable compression injury.

Conclusion: In the edematous tissue condition, unintentional overcompression could be increased and result in tissue injury on the compression line of the circular stapler.

[Ann Surg Treat Res 2020;99(2):72-81]

Key Words: Anastomotic leak, Colorectal surgery, Soft tissue injuries, Surgical anastomosis, Surgical staplers

INTRODUCTION

The proportion of sphincter-preserving surgery for rectal

cancer has increased dramatically in the past several decades with the development of the 'double stapling technique' since both circular and linear staplers allow for colorectal continuity

Received March 10, 2020, Revised May 12, 2020, Accepted May 15, 2020

Corresponding Author: Gyung Mo Son

Department of Surgery, Pusan National University Yangsan Hospital, 20 Geumo-ro, Mulgeum-eup, Yangsan 50612, Korea

Tel: +82-55-360-2124, **Fax:** +82-55-360-2154

E-mail: skm1711@pusan.ac.kr

ORCID: https://orcid.org/0000-0002-8861-6293

Copyright © 2020, the Korean Surgical Society

© Annals of Surgical Treatment and Research is an Open Access Journal. All articles are distributed under the terms of the Creative Commons Attribution Non-Commercial License (<http://creativecommons.org/licenses/by-nc/4.0/>) which permits unrestricted non-commercial use, distribution, and reproduction in any medium, provided the original work is properly cited.

in the deep pelvic cavity [1,2]. The advancement of stapling devices has enabled safe standardized stapled anastomosis in clinical practice [3,4]. The goal of good stapling technique is to achieve tissue apposition and hemostasis without tissue damage and ischemia [5].

Despite the technical improvements, adverse events by circular stapler, including tissue damage, staple line defect, and malformed staples, have been reported in 19% in the clinical field [5]. Particularly, tissue damage such as seromuscular rupture is not uncommon in the compression process when the colon wall had been weakened by colon cancer obstruction or preoperative chemoradiation therapy for rectal cancer [2,6]. Recently, the U.S. Food and Drug Administration announced that over 41,000 individual medical device reports for surgical staplers accumulated during the last 7 years including 366 deaths, over 9,000 serious injuries, and over 32,000 malfunctions. As such, they have emphasized the importance of the proper use of surgical staplers to prevent surgical stapler-related casualties and serious malfunctions and improve patient safety [7].

Tissue compression is an important step to ensure safe staple formation [8,9]. The degree of proper compression is affected by variable tissue conditions such as tissue thickness, density, location, edema, and pathologic status [8]. Surgeons usually consider these kinds of tissue conditions to obtain an optimal compression. However, such efforts sometimes do not completely eliminate the chance of accidental tissue damage during the compression process by circular stapler.

Improper usage of a circular stapler may result in tissue injury that could be considered an important risk factor for anastomotic leak [3,8]. However, it is not easy for a surgeon to find tissue damage in the standardized stapling process. Especially, a crushing tissue injury on the compression line is often overlooked by surgeons because the stapled colon wall has been inverted to the luminal side which causes difficulty for the operator to notice the stapled tissue damage [2]. That is why such tissue damage is not often identified during the rectal cancer surgery in the deep pelvic cavity. Thus, surgeons have not been interested in tissue compression injury due to circular staples, until recently. In addition, due to the lack of studies related to tissue compression injury due to circular staples, its damaging mechanism was not clearly identified. Therefore, the authors designed a preliminary *in-vitro* study to quantitatively evaluate compression injury of circular stapling line on artificial collagen tissue.

This study was aimed to evaluate the risk factors and explore the mechanism of tissue compression injury from the use of a circular stapler for end-to-end anastomosis using *in-vitro* model.

METHODS

Preparation for experiments

This study was conducted with the approval of the Institutional Animal Care and Use Committee of Pusan National University Yangsan Hospital (No. H-2018-060). For the visual assessment of compression injury *in-vitro* model, transparent collagen plates were produced by mixing 80% of porcine collagen with 20% of glycerin in a rectangular shape with a thickness of 2 and 3 mm (collagen plate; SINI Inc., Uiwang, Korea). The edematous collagen plates were prepared by tap water immersion for 2 minutes.

Physical property measurements

To measure the physical properties, both dry and edematous wet collagen plates were compared with fresh colon tissues from healthy porcine that were slaughtered within 24 hours. The tensile strength (g/mm^2) was measured at a table speed 120 mm/min and load cell with a maximum of 10 kg in weight. The deformation (mm) was measured at a table speed of 60 mm/min and load cell with a maximum of 2 kg in weight. The springiness (%) and cohesiveness (%) were measured by a double compression test of tissue texture. The physical properties were examined by using rheometer (Type COMPAC-100II; Sun Science Co., Tokyo, Japan).

In-vitro experiments

Two commercial circular stapler devices were used in the *in-vitro* experiments: (1) adjustable-compression type stapler: CDH29A, 29-mm diameter (Ethicon Endo-Surgery, Inc., Cincinnati, OH, USA) with adjustable closed staple height (1–2.5 mm); and (2) fixed-compression type stapler: EEA, 28-mm diameter (Covidien, Inc., Mansfield, MA, USA) with green cartridge at 2 mm of closed staple height. Once the firing trigger was fully closed, the heights of the staple driver of the CDH and pusher bar of the EEA as a homologous unit were measured using a digital caliper (Digimatic caliper, No.500-151-30; Mitutoyo, Kanagawa, Japan). The shape of the anvil was also compared between the 2 circular staplers (Fig. 1A, B). The dry and wet collagen plates were prepared to reproduce normal and edematous tissue conditions. Two collagen plates with the same thickness and tissue conditions were placed on the staple cartridge through the anvil needle (Fig. 1C). The examiners rotated the compression knob with 2 fingers for anvil approximation toward the staple housing.

Compression type

In the adjustable-type stapler, controlled tissue compression (CTC) scale on the body shaft was used to measure the gap distance as closed staple height ranging 1.0–2.5 mm (Fig. 1C, D). For objective evaluation, CTC scale was covered with opaque

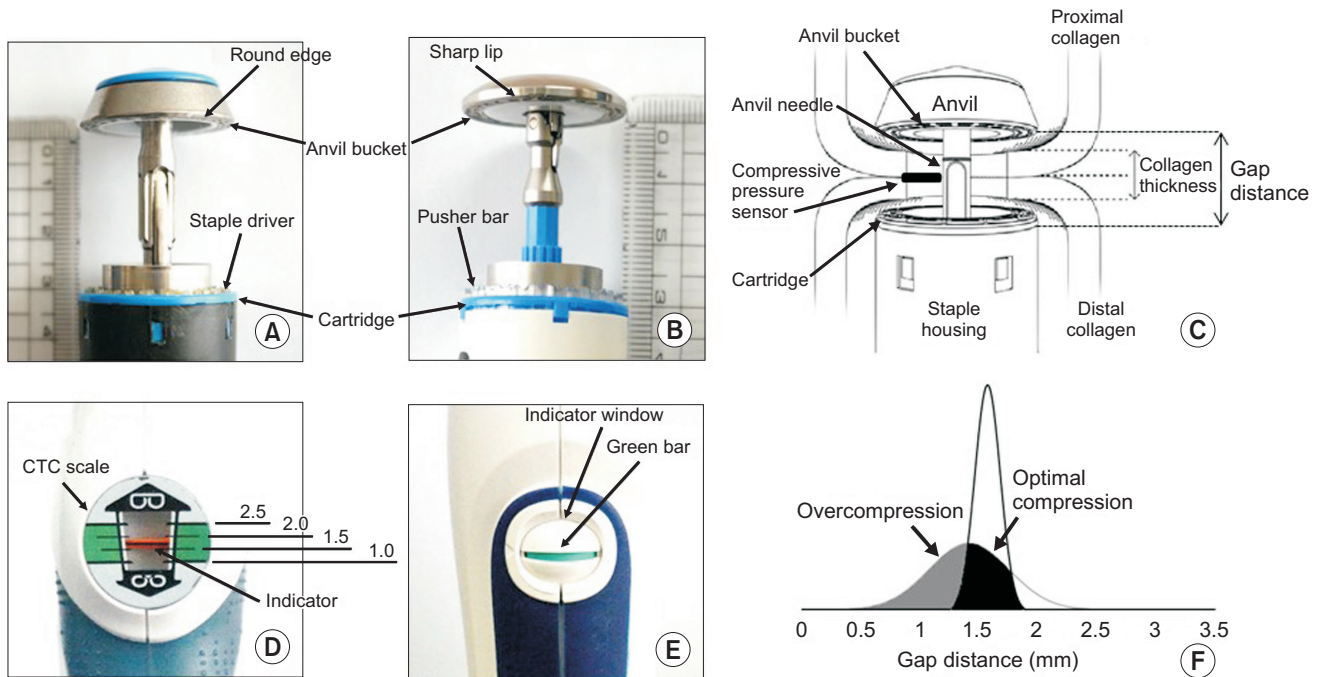


Fig. 1. Structure of 2 commercial circular stapler devices. (A) Anvil connecting unit of adjustable-compression type stapler showed round anvil edge and small staple driver on the staple cartridge when firing trigger was fully closed. (B) Fixed-compression type stapler had sharp lip on anvil edge and pusher bar protruding from the staple cartridge. (C) The cross-sectional view of collagen plates on the circular stapler with compressive pressure sensor placed between proximal and distal collagen plates. (D) Controlled tissue compression (CTC) scale and (E) indicator window with green bar on the stapler shafts using indicators for proper compression. (F) Optimal and overcompression depending on the gap distance between anvil and staple cartridge.

paper and the thickness and tissue condition of the collagen plates were blinded to the examiners. To perform 'adjustable compression,' the examiners used tactile feedback of sudden increase of resistance between the 2 fingers. Once compression was applied for 15 seconds, level of gap distance on CTC scale was recorded. Additionally, the examiner applied the 'maximal compression' by turning the rotating knob to the minimum gap distance (1 mm). In the fixed-type stapler, the anvil was approximated until the green bar was completely visualized on the indicator window of stapler's shaft, and it was categorized as 'fixed compression' (Fig. 1E). The gap distance of the fixed-type stapler was measured between anvil edge and staple cartridge through transparent collagen plates using a digital caliper. In the initial pilot experiment of this study, the dry collagen plates were compressed without serious crushing injury. Therefore, 'optimal compression' was defined as compression within the gap distance range of dry collagen. If the gap distance was shorter than the optimal compression range, it was considered as 'overcompression' (Fig. 1F).

Compressive pressure measurement

The real-time compressive pressure was measured with a commercial pressure sensor (FlexiForce B201-H; Tekscan, Inc., Boston, MA, USA). The compressive pressure sensor was

inserted between the proximal and distal collagen plates and compressive force was measured in the anvil approximation process (Fig. 1C). The compressive pressure was digitalized to load (N) using 200 Hz of sampling rate and displayed by ELF System ver. A software (Tekscan, Inc).

Compression injury scale

The tissue was mounted on a circular stapler and compressed for 15 seconds. The compressed collagen plates were removed carefully from the anvil needle of the circular stapler and compression injury on the proximal and distal collagen plates was accessed. The depth of injury was evaluated through a cross section of the transparent collagen plates by 2 independent examiners under blinded states of collagen condition and compression type. The compression injury scale was designed by the authors based on the Likert scale [10]. Acceptable injury was defined as cases in which the depth of tissue injury extended to less than half of collagen thickness; unacceptable injury was defined as cases where the damage exceeded half of the collagen thickness (Fig. 2).

Statistical analysis

The tissue factors were included as collagen thickness (4 or 6 mm), tissue condition (dry or wet), and side location (proximal

or distal). The compression factors were compression type (adjustable, maximal, or fixed) and compression propriety (optimal or overcompression). The unpaired t-test was used to compare the mean value of continuous variables. Chi-square test or Pearson correlation test was applied to compare the categorical variables. Physical property data of collagen plates and porcine colon tissue were compared using 1-way analysis of variance (ANOVA) and the Bonferroni test was used for *post-hoc* analysis. Bivariate analysis was used for correlation analysis between significant risk factors. Multivariate analysis

was performed with binary logistic regression model using a forward condition analysis. The covariance input criterion was less than 0.1 and elimination criterion was less than 0.05. Statistical analysis was performed using IBM SPSS Statistics ver. 24.0 (IBM Corp., Armonk, NY, USA) and $P < 0.05$ was considered as significance level.

RESULTS

The physical properties of the collagen plates were compared with those of porcine colon tissue (Table 1). The tensile strength and tissue texture such as deformation, springiness, and cohesiveness were comparable between the dry collagen plates and colon tissue. In the wet collagen plates, the tensile strength was significantly weaker than dry collagen or healthy colon tissue.

There were some structural differences between the 2 types of staplers. The height of pusher bar in the fixed-type stapler was significantly higher than those of the staple driver in the adjusting-type stapler (1.95 ± 0.01 mm vs. 0.73 ± 0.03 mm, $P < 0.001$). In the shape of anvil edge, a sharp lip protruded around the anvil of the fixed-type stapler and the length of sharp lip was 0.30 ± 0.01 mm. The adjustable-type stapler had an atraumatic round shape of the anvil edge.

The outcomes of *in-vitro* experiments were different according to the collagen thickness and tissue conditions (Table 2). In the adjustable compression, the compressibility ratio was similar between dry and wet collagen plates (61.7% vs. 63.2%, respectively, $P = 0.234$). However, on the edematous collagen plates, the standard deviation ranges of gap distance were spread out more than 3 times those of dry plates. The rate of overcompression and unacceptable injury also had increased in the wet collagen. The overcompression cases had significantly lower gap distance and higher incidence of

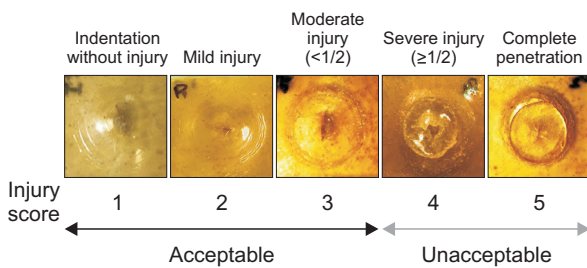


Fig. 2. Compression injury scale. Acceptable injury was defined as cases in which the depth of tissue injury extended to less than half of collagen thickness (score 1–3); unacceptable injury was defined as cases where the depth of damage exceeded the half of collagen thickness (score 4–5).

Table 1. Physical properties of dry and wet collagen plates compared with fresh porcine colon tissue

Variable	Dry collagen	Wet collagen	Porcine colon
Tensile strength (g/mm ²)	385.0 ± 7.1	184.0 ± 5.7*	387.5 ± 10.6
Deformation (mm)	2.0 ± 0.0	2.0 ± 0.1	3.2 ± 0.1
Springiness (%)	105.3 ± 3.0	92.3 ± 7.3	93.6 ± 4.8
Cohesiveness (%)	96.5 ± 3.6	92.9 ± 0.9	92.3 ± 7.3

* $P < 0.05$.

Table 2. The gap distance, overcompression, and unacceptable injury rate according to the compression type

Compression type	Collagen plates			Gap distance	Overcompression	Unacceptable injury
	Thickness (mm)	Condition	Number			
Adjustable compression	4	Dry	48	1.53 ± 0.08 (1.4–1.6)	-	2 (4.2)
		Wet	68	1.47 ± 0.36 (1.1–2.0)	32 (47.1)	9 (13.2)
	6	Dry	48	2.35 ± 0.13 (2.1–2.5)	-	0 (0)
		Wet	68	2.22 ± 0.38 (1.6–2.5)	24 (35.3)	8 (11.8)
Maximal compression	4	Dry	48	1.01 ± 0.03 (1.0–1.1)	-	0 (0)
		Wet	68	1.03 ± 0.05 (1.0–1.1)	68 (100)	9 (13.2)
	6	Dry	48	1.02 ± 0.04 (1.0–1.1)	-	1 (2.1)
		Wet	68	1.04 ± 0.05 (1.0–1.1)	68 (100)	22 (32.4)
Fixed compression	4	Dry	20	1.03 ± 0.18 (0.9–1.3)	-	0 (0)
		Wet	24	0.68 ± 0.11 (0.5–0.8)	24 (100)	15 (62.5)
	6	Dry	20	1.52 ± 0.19 (1.2–1.8)	-	1 (5.0)
		Wet	24	0.95 ± 0.09 (0.8–1.1)	24 (100)	13 (54.2)

Values are presented as mean ± standard deviation (range) or number (%).

unacceptable injuries than those of optimal compression cases in the wet collagen plates (Fig. 3). By maximal compression of the adjustable stapler, gap distance was intentionally made to 1

mm. Unacceptable injury was little on the dry collagen plates, but it increased to 32.4% on the wet, 6-mm thick collagen plates. In the fixed compression, gap distance was significantly

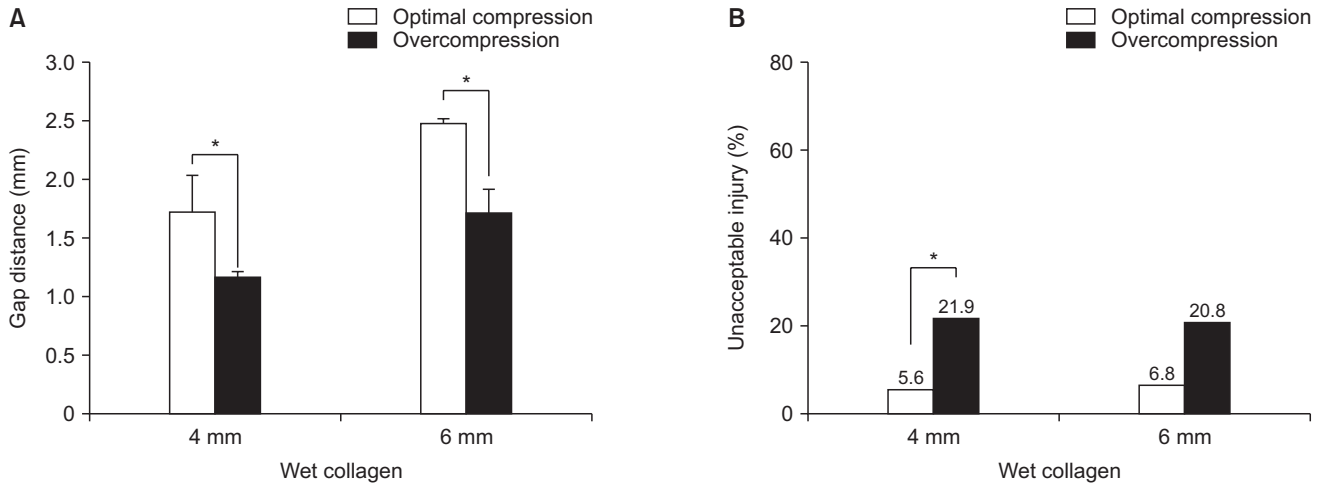


Fig. 3. The gap distance and unacceptable injury of the adjustable type stapler according to compression propriety on the wet collagen plates. (A) Comparing with optimal compression, the overcompression cases had significantly lower gap distance. (B) The incidence of unacceptable injuries was higher in the overcompression cases than optimal compression cases. * $P < 0.05$.

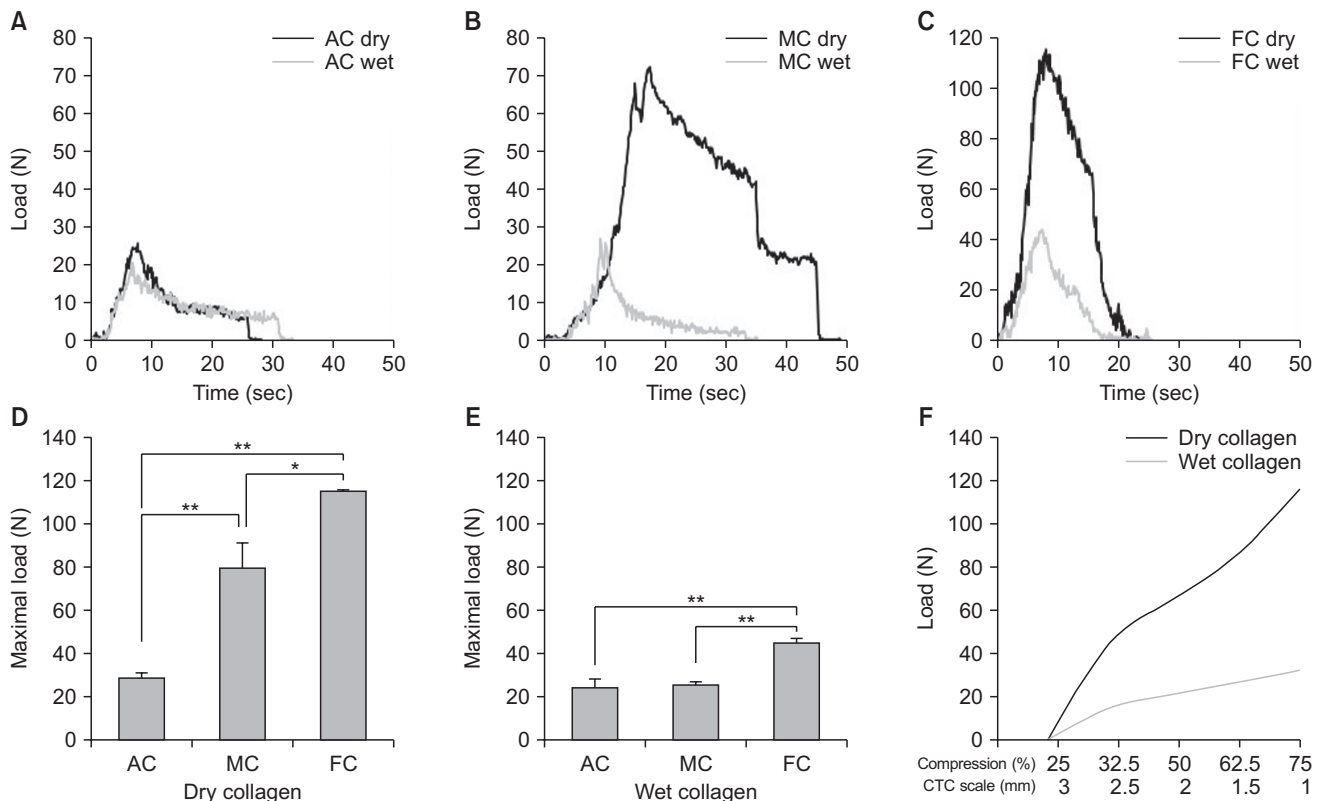


Fig. 4. Real-time compressive pressure measurements. Compressive pressure curves according to compression type of AC (A), MC (B), and FC (C) on the both dry and wet collagen plates. Peak compressive pressure levels showed significant increment in FC on the dry (D) and wet (E) collagen plates. (F) The upward curve patterns of compressive pressure were different on the collagen conditions. In contrary to steep pressure rising on the dry collagen plates, gradual plateau was showed even after optimal compression point on the wet collagen. AC, adjustable compression; MC, maximal compression; FC, fixed compression; CTC scale, controlled tissue compression scale. * $P < 0.05$, ** $P < 0.001$.

narrower than that of adjustable compression ($P < 0.001$). The compressibility ratio increased significantly from dry collagen to wet condition (74.2% vs. 83.1%, respectively, $P < 0.001$). Unacceptable injuries were observed in more than 50% of the wet collagens cases whereas only compressive impression or slight injuries were found on the dry collagen plates (Table 2).

In the real-time compressive pressure measurements, the adjustable compression showed similar pressure patterns between the dry and wet collagen plates (Fig. 4A). In the maximal or fixed compression, compressive pressures were markedly increased in the dry collagen plates than in the edematous collagen plates (Fig. 4B, C). Peak compressive pressure levels showed significant increments in the fixed compression than those of adjustable compression in both dry and wet collagen plates (Fig. 4D, E). As the gap distance was approximated in the adjustable compression, the upward curve patterns of compressive pressure were different depending on the collagen conditions. On the dry collagen plates, the compressive pressure curve showed a steep increase as the gap distance was closing. In contrast, the compressive pressure reached a gradual plateau even after reaching an optimal compression point on the wet collagen plates (Fig. 4F).

On bivariate analysis of the risk factors relating to unacceptable injury, the fixed compression and wet conditioned plates were respectively correlated with overcompression which was considered as the prerequisite causing unacceptable injury (Fig. 5).

On univariate analysis, the rates of unacceptable injury were significantly higher in the wet collagen plates, proximal anvil side, fixed compression, and overcompression. On multivariate analysis, edematous collagen condition was the most important risk factor. The proximal anvil side, fixed compression, and overcompression were also analyzed as independent risk factors of unacceptable injury (Table 3).

DISCUSSION

Despite the development of stapler devices, it is not uncommon that tissue injury occurs on the staple line of edematous or friable bowel in the surgical field [3,11,12]. In this study, both adjustable and fixed compression were performed safely without damage on the dry collagen plates as normal healthy tissue. However, on the edematous wet plates, the rates of unintentional overcompression were up

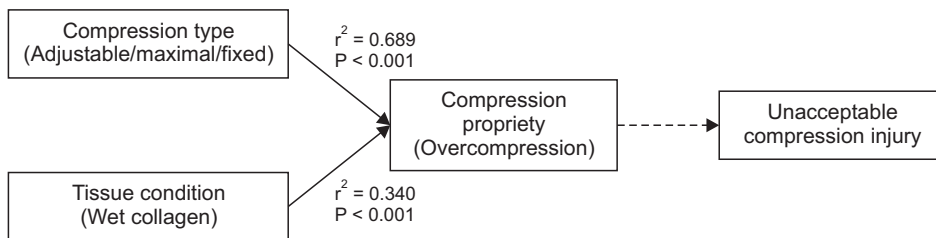


Fig. 5. The bivariate correlation analysis of the risk factors relating unacceptable injury.

Table 3. Univariate and multivariate analysis for the risk factors of compression injury

Risk factor	Number	Unacceptable injury		Multivariate analysis		
		No. (%)	P-value	Odds ratio	95% CI	P-value
Collagen condition						
Dry	232	4 (1.7)	<0.001	1.000		
Wet	320	76 (23.8)		13.919	4.535–42.719	<0.001
Thickness						
4 mm	276	35 (12.7)	0.227			
6 mm	276	45 (16.3)				
Side						
Distal	276	22 (8.0)	<0.001	1		
Proximal	276	58 (21.0)		4.234	2.314–7.746	<0.001
Compression type						
Adjustable	232	19 (8.2)	<0.001	1		
Maximal	232	32 (13.8)		0.997	0.463–2.151	0.995
Fixed	88	29 (33.0)		5.794	2.494–13.461	<0.001
Compression propriety						
Optimal	216	8 (3.7)	<0.001	1		
Overcompression	336	72 (21.4)		3.671	1.394–9.669	0.008

CI, confidence interval.

to 35%, and approximately 30% of them were associated with serious compression injury from adjustable compression. Fixed compression also caused overcompression in all cases, and tissue damage was observed in more than 50% of wet collagen plates. Even though there were different outcomes for tissue damage, edematous condition could be considered as the most significant risk factor for compression injury rather than the compression types of circular stapler device. Therefore, compression injury might be unavoidable on edematous tissue. However, understanding its mechanism of compression injury may help surgeons find the proper and safe usage of circular staplers.

In this study, the wet collagen plates and compression types were significantly related with overcompression which is considered an important prerequisite causing unacceptable compression injury. In order to evaluate the mechanism of the overcompression phenomenon, we measured the compression pressure patterns in proportion to the approximation of gap distance. In the adjustable compression device, the initial compression pressure increased until the gap distances reached almost 2 mm. The compression pressures vary between 60–70 N and these values were close to 6–8 g/mm², considering the target compressed area was 10 mm² in this model. Such compressive force is well associated with the stapler manufacturer's design concept to deliver 6–8 g/mm² of compressive force to the target area based on the previous researches seeking optimal compression [8,13,14].

In the dry collagen plates, a steep increase of compressive pressure was observed as gap distance reached the optimal point. This sudden change of compressive pressure serves as a useful indicator for the operator to recognize its optimal compression level. Healthy tissue has both a structural stability and viscoelasticity, which reflects its physical properties of high density with good elasticity. As the tissue is compressed slowly, the resistance force acting on the tissue surface against compression force is also slowly increased since the intercellular space inside the tissue is shrunken [3]. This is a biomechanics perspective of how force is acting on tissue upon compression. Once the healthy tissue is gently, but completely, compressed, when the intercellular fluid is spilled out and tissue stiffness increase without cellular damage is established, the resistance force reaches the strongest point. The operator can perceive the abrupt resistance triggering turning of the rotating knob of the circular stapler device.

However, in the edematous collagen plates, the compressive pressure curve formed a gradual plateau even after further compression was applied beyond the optimal compression point. If there were presence of edema or inflammation in the colon, the tissue structure could be weakened and its density, elasticity, and restoration capability would have been deteriorated [15]. As a result, its structural integrity could be

even further weakened by the accumulation of muscular and submucosal tears caused by overcompression [12,16]. In the edematous tissue, its internal structure is not as stable as that of normal tissue and sufficient resistance force could not establish action on the tissue surface against compression force [17,18]. In this case, compressive pressure curve formed a plateau curve passing through the optimal point, even if compression is continued.

Thus, even an experienced operator could not often judge the optimal compression point on the edematous tissue using tactile feedback of compressive pressure, and it is easy to pass the optimal compression point inducing overcompression [17,18]. Therefore, the range of standard deviation for the gap distance was significantly wider in the wet collagen plates than that of dry collagen in this experiment. If excessive compression is applied beyond the optimal point, tissue damage could occur more frequently. That is why we have proposed this overcompression phenomenon as one of the key potential mechanisms causing compression injury in the edematous tissue. If surgeons fail to recognize the optimal compression point and apply overcompression, tissue damage could increase the risk of compression tissue injury and anastomotic complications [3]. Therefore, a surgeon should be fully familiar with the proper usage of the circular stapler according to different tissue conditions since inappropriate usage could lead to a serious threat to patient safety [3,17]. Also, we want to suggest that current adjustable compression based on the resistance feedback system may be suboptimal in preventing overcompression in edematous and friable tissue conditions.

The fixed compression device is designed to present a green bar indicator window of proper compression. Operators usually rotate the knob until the green bar is visualized, and that is the only indicator that proper compression has been achieved and the device is ready to fire. The gap distance of the fixed-compression type device was measured as shorter, at 0.5–0.8 mm, than the adjustable-type device on the dry collagen plates. There is a different method of measuring gap distance; using the CTC scale for the adjustable-type stapler and digital caliper for the fixed-type stapler. As such, almost 0.5 mm of anvil bucket and 0.3 mm of sharp lip on the anvil edge unit should be considered as an additional gap that might form similar gap distance ranges to the adjustable-type device. This could explain the comparable tissue injury rates on the dry collagen plates with the different compressibility of the 2 kinds of circular staplers.

Interestingly, the fixed-type compression caused more frequent tissue damage on the edematous tissue. In this study, researchers looked for additional causes of tissue damage from the structural differences between the fixed-type and adjustable-type devices. These different structural features of the circular staplers might have caused more frequent tissue

damage on the edematous tissue. First, the fixed-type device had almost 10% excessive compression on the wet collagen plates than dry plates, but the adjustable-type device showed a similar level of compressibility between dry and wet collagen plates. Second, the fixed compression device has a sharp lip on the anvil edge while the adjustable device has an atraumatic round-shaped anvil. This edged structure appears to form dissipative frictional contact to limit tissue slippage in a radial direction during anvil approximation, which could prevent tissue elongation maintaining original tissue thickness and might have function in limiting maximal tissue compression [18]. Unfortunately, this anvil's sharp lip could cause shearing on the edematous tissue. Third, there were different sized staple stationary units in the cartridges of the fixed and adjustable-type devices. The cartridges of fixed and adjustable-type instruments are equipped with a staple driver or pusher bar as homologous devices for staple stationary units. This staple driver holds the staple inside the cartridge and pushes staples towards the anvil bucket to form a B shape when the firing trigger is closed. However, in the fixed-type device, the pusher bar was more elevated than the staple driver of the adjustable-type devices when the trigger was fully closed. These different structural features of circular staplers might have caused more frequent tissue damage on the edematous tissue.

On the dry collagen conditions which represent normal tissue, the performance of both circular stapler devices was safe and acceptable. However, with edematous or inflammatory conditions, the different device structures could act as a dangerous trap with improper usage in the compression process. Therefore, the authors would like to suggest that structural aspects relating to the potential risk of tissue damage should be improved upon in the development of new circular staples.

This study had several limitations. First, the single-layered collagen plate is a simple structure compared to multi-layered colon tissue. It is possible that compression injury in the collagen plate may be exaggerated when it is compared to edematous colon tissue [19]. On the physical property test, tensile strength was similar between dry collagen and healthy colon tissue. Thus, unacceptable compression injuries were rarely seen in the dry collagen plate regardless of compression type. However, the edematous collagen plate had significantly lower tensile strength than those of the dry collagen and healthy colon tissue. Therefore, the wet collagen plate was considered to be more susceptible to crushing injury. However, the collagen model differs from actual colon tissue in several aspects. The collagen model has a single layer structure of mixed components of collagen and glycerin, whereas the large intestine tissue has a multilayer structure with different thicknesses and physical properties. Therefore, the degree of tissue compression and its damage on the collagen plates

must also be different from the actual colon tissue. Therefore, the results of this *in-vitro* study need to be verified in real colon tissues. In addition, since animal colon and human colon tissues also have different thicknesses and strengths, subsequent studies need to be continued using human colon tissues obtained during surgery. In spite of these structural differences, collagen plates have some advantages on simulating edematous tissue conditions simply by soaking it in tap water for 2 minutes. Also, the transparent plate makes it possible to provide visual feedback of tissue injury and to assess the degree of crushing damage. This collagen plate model is also inexpensive, easy to manufacture, and free from ethical constraints which make it possible to carry out repetitive exercises for educational information regarding the safe use of circular staplers. Therefore, the collagen plates could be considered an effective training tool in learning how to determine optimal compression of circular staplers for novices using a circular stapler in real practice. Second, this study did not assess the tissue perfusion and bursting strength of stapled tissue. Optimal tissue compression should avoid tissue damage and maintain adequate tensile strength with adequate blood flow [20]. Third, we performed only full compression of the fixed-type stapler until the green bar became fully visualized in this study. Therefore, experiments using partial visualization as adjustable compression were not conducted for EEA devices in this study. The instability of fixed-type staplers may have originated from specifically determinant ways of using complete visualization of the green bar for proper compression. The instructions for users (IFU) of EEA circular staplers described that the green bar must be visible in the indicator window before releasing the safety lever prior to firing. The in-service educational information of EEA device recommends that full visualization of the green bar indicates that the stapler is ready to be fired; which is generally used as the proper compression indicator during surgery. However, even if only a part of the green bar is visible, the EEA device is designed to enable firing. The authors had some experiences of safe stapling of EEA when the green bar was partially visible on thick colon tissues in the surgical field. This partial visualization related patent was registered and that novel indicator window was designed to enable adjustable compression, but this patent technology has not been applied to the commercial EEA device until now [21]. Therefore, even if the fixed-type stapler is used, firing may be attempted when colon tissue can be compressed comfortably at appropriate compression levels. However, since the indicator window for viewing the compression scale is narrow and there is no detailed mention of using partial or full visualization of the green bar on the indicator window for optimal compression in the IFU of the EEA device, inexperienced users still could have uncertainty in achieving optimal compression using circular staplers. Therefore, the researchers eagerly hope for

the advancement and development of circular stapler devices in which the operator can easily and clearly recognize optimal compression levels using a distinctive indicator.

This was a preliminary *in-vitro* study on the effect of compression injury of the circular stapling line on artificial collagen tissue. Therefore, the results based on quantitative evaluation of compression injury suggest additional *ex-vivo* evaluations using real colon tissue in the future. Clinical studies also should be conducted to identify the clinical implications of crushing injuries caused by overcompression of circular staplers. Also, further clinical studies should be designed to help develop effective training programs for the standardization and optimization for proper application of staple devices.

In edematous tissue conditions, unintentional overcompression could be increased and result in unacceptable tissue injury on the compression line of circular staplers. Well-designed training programs for new device users are required to teach proper usage of circular staplers to minimize compression injury at the stapling line to ensure safe anastomosis.

ACKNOWLEDGEMENTS

The authors thank to Susung Kim and Min Jung Choi for technical assistance and Michael Mun Ki Choi for illustration.

Fund/Grant Support

This study was supported by Research Institute for

Convergence of Biomedical Science and Technology Grant (30-2018-025), Pusan National University Yangsan Hospital.

Conflict of Interest

No potential conflict of interest relevant to this article was reported.

ORCID iD

Gyung Mo Son: <https://orcid.org/0000-0002-8861-6293>

Myeong Sook Kwon: <https://orcid.org/0000-0003-3433-6744>

Hong-min Ahn: <https://orcid.org/0000-0001-9963-2021>

In Young Lee: <https://orcid.org/0000-0001-5954-6188>

Gun Ho Kim: <https://orcid.org/0000-0001-8281-1701>

Kyoung Won Nam: <https://orcid.org/0000-0002-8378-8758>

JoonWon Lee: <https://orcid.org/0000-0002-4102-7440>

Chang-Suk Kong: <https://orcid.org/0000-0002-2304-9304>

Author Contribution

Conceptualization: GMS

Formal Analysis: GMS, MSK, IYL

Investigation: GMS, MSK, HA

Methodology: GHK, KWN, JWJ, CSK

Project Administration: GMS, MSK

Writing – Original Draft: GMS, MSK, GHK

Writing – Review & Editing: GMS, HA, IYL, KWN, JWJ, CSK

REFERENCES

- Nakayama S, Hasegawa S, Hida K, Kawada K, Sakai Y. Obtaining secure stapling of a double stapling anastomosis. *J Surg Res* 2015;193:652-7.
- Ikeda T, Kumashiro R, Taketani K, Ando K, Kimura Y, Saeki H, et al. Endoscopic evaluation of clinical colorectal anastomotic leakage. *J Surg Res* 2015;193:126-34.
- Offodile AC 2nd, Feingold DL, Nasar A, Whelan RL, Arnell TD. High incidence of technical errors involving the EEA circular stapler: a single institution experience. *J Am Coll Surg* 2010;210:331-5.
- Slieker JC, Daams F, Mulder IM, Jeekel J, Lange JF. Systematic review of the technique of colorectal anastomosis. *JAMA Surg* 2013;148:190-201.
- Mery CM, Shafi BM, Binyamin G, Morton JM, Gertner M. Profiling surgical staplers: effect of staple height, buttress, and overlap on staple line failure. *Surg Obes Relat Dis* 2008;4:416-22.
- Rho SY, Bae SU, Baek SJ, Hur H, Min BS, Baik SH, et al. Feasibility and safety of laparoscopic resection following stent insertion for obstructing left-sided colon cancer. *J Korean Surg Soc* 2013;85:290-5.
- Maisel W. Safe use of surgical staplers and staples - letter to health care providers [Internet]. Silver Spring: US Food and Drug Administration; 2019 [cited 2020 Jul 3]. Available from: <https://www.fda.gov/medical-devices/letters-health-care-providers/safe-use-surgical-staplers-and-staples-letter-health-care-providers>.
- Chekan E, Whelan RL. Surgical stapling device-tissue interactions: what surgeons need to know to improve patient outcomes. *Med Devices (Auckl)* 2014;7:305-18.
- Nakayama S, Hasegawa S, Nagayama S, Kato S, Hida K, Tanaka E, et al. The importance of precompression time for secure stapling with a linear stapler. *Surg Endosc* 2011;25:2382-6.
- Sureshkumar S, Jubel K, Ali MS, Vijayakumar C, Amaranathan A, Sundaramoorthy S, et al. Comparing surgical site infection and scar cosmesis between conventional linear skin closure versus purse-string skin closure in stoma reversal - a randomized controlled trial. *Cureus* 2018;10:e2181.
- Matsuzawa F, Homma S, Yoshida T.

- Konishi Y, Shibasaki S, Ishikawa T, et al. Serosal laceration during firing of powered linear stapler is a predictor of staple malformation. *Surg Innov* 2017;24:590-7.
12. Ikeda T, Kumashiro R, Oki E, Taketani K, Ando K, Aishima S, et al. Evaluation of techniques to prevent colorectal anastomotic leakage. *J Surg Res* 2015;194:450-7.
13. Baker RS, Foote J, Kemmeter P, Brady R, Vroegop T, Serveld M. The science of stapling and leaks. *Obes Surg* 2004;14:1290-8.
14. Yang J, Soltz M, Russell H, Beres J, Zhao J, Liao D, et al. Surface deformation analysis of end-to-end stapled intestinal anastomosis. *Surg Innov* 2012;19:281-7.
15. Kimura M, Takahashi H, Tanaka T, Sagawa H, Takeyama H. Profiling surgical staplers: does staple direction affect the strength of the anastomosis. *Dig Surg* 2015;32:454-8.
16. Egorov VI, Schastlivtsev IV, Prut EV, Baranov AO, Turusov RA. Mechanical properties of the human gastrointestinal tract. *J Biomech* 2002;35:1417-25.
17. Myers SR, Rothermel WS Jr, Shaffer L. The effect of tissue compression on circular stapler line failure. *Surg Endosc* 2011;25:3043-9.
18. McGuire J, Wright IC, Leverment JN. An in vitro assessment of tissue compression damage during circular stapler approximation tests, measuring expulsion of intracellular fluid and force. *Proc Inst Mech Eng H* 2001;215:589-97.
19. Watters DA, Smith AN, Eastwood MA, Anderson KC, Elton RA, Mugerwa JW. Mechanical properties of the colon: comparison of the features of the African and European colon in vitro. *Gut* 1985;26:384-92.
20. Chung RS. Blood flow in colonic anastomoses: effect of stapling and suturing. *Ann Surg* 1987;206:335-9.
21. Zhang X, inventor; Covidien LP, assignee. Circular stapler with visual indicator mechanism. United States patent US 16/073,354. 2019 Feb 7.

©2020 Ethicon, Inc. All rights reserved. 149606-200812

Ethicon Inc. has no independent knowledge concerning the information contained in this article, and findings and conclusions expressed are those reached independently by the authors