

HEALIX ADVANCE™

Knotless Anchors

The HEALIX Family of Anchors™ provides a comprehensive and versatile system for rotator cuff repairs.



HEALIX ADVANCE™ Anchors Technology¹

- Multi-thread design engages cortical and cancellous bone effectively¹
- Secure fixation demonstrated through high pullout strength²
- Cannulated to channel blood flow to the surface

Versatility

Versatile Suture Configurations

Size & material

- Implants available in 4.75, 5.5, and 6.5 mm sizes.
- Implants available in PEEK and BIOCRYL RAPIDE® Biocomposite Material.

Dual Wire Threader

Allows increased suture passes through the anchor

- Loads up to 6 strands of #2 ORTHOCORD® Suture/ PERMACORD™ Suture (4 proximal 2 distal)
- 5 strands of #2 DYNACORD™ Suture
- 3 strands of PERMATAPE™ Suture

Usability

The Sliding Suture Ring

Suture wrapping

- Prevents suture wrapping inside the cannula, while automatically unwinding stitches from above the ring as the inserter is removed.

Suture tensioning

- 360° cleats enable suture management during tensioning.

Stay-Stitch Sliding #2 ORTHOCORD® Suture

- Can be incorporated into the repair.

Locking suture ring


- For open use



Product Codes

HEALIX ADVANCE Anchors

Item Code	Item Description
222885	4.75mm HEALIX ADVANCE Knotless BR Anchor
222886	5.5mm HEALIX ADVANCE Knotless BR Anchor
222887	6.5mm HEALIX ADVANCE Knotless BR Anchor
222888	4.75mm HEALIX ADVANCE Knotless PEEK Anchor
222889	5.5mm HEALIX ADVANCE Knotless PEEK Anchor
222890	6.5mm HEALIX ADVANCE Knotless PEEK Anchor

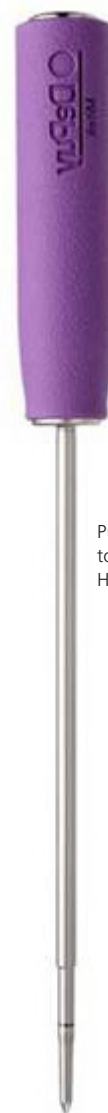

Purple handle (222314), to be used with 4.75 mm HEALIX ADVANCE Knotless Anchors.

HEALIX ADVANCE Instrumentation

Item Code	Item Description
222314	4.5 HEALIX ADVANCE Anchor and 4.75 HEALIX ADVANCE Knotless Anchor Awl
222315	5.5/6.5 HEALIX ADVANCE Anchor Awl
222223	5.5/6.5 HEALIX ADVANCE Knotless Anchor Awl
222316	4.5 HEALIX ADVANCE Anchor Cortical Awl/Tap
222317	5.5/6.5 HEALIX ADVANCE Anchor Cortical Awl/Tap
222328	4.5 HEALIX ADVANCE Anchor Full Tap
222329	5.5/6.5 HEALIX ADVANCE Anchor Full Tap



Grey handle (222223), to be used with 5.5 and 6.5 mm HEALIX ADVANCE Knotless Anchors.

References:

1. Internal test report RR-101380-R Rev A ACS
2. Barber FA, Herbert MA, Beavis RC, Barrera Oro FB. Suture anchor materials, eyelets and designs. update 2008. Arthroscopy. 2008;24(8):859.867.

Limited Warranty and Disclaimer: DePuy Synthes Mitek Sports Medicine products are sold with a limited warranty to the original purchaser against defects in workmanship and materials. Any other express or implied warranties, including warranties of merchantability or fitness, are hereby disclaimed.

WARNING: Please also refer to the package insert(s) or other labeling associated with the devices identified in this brochure for additional information.

This publication is not intended for distribution in the USA.

The third party trademarks used herein are the trademarks of their respective owners.

For full product details and precautions please consult the IFU.

