


HEALIX ADVANCE™

Knotless Anchors

The HEALIX Family of Anchors™ provides a comprehensive and versatile system for rotator cuff repairs.

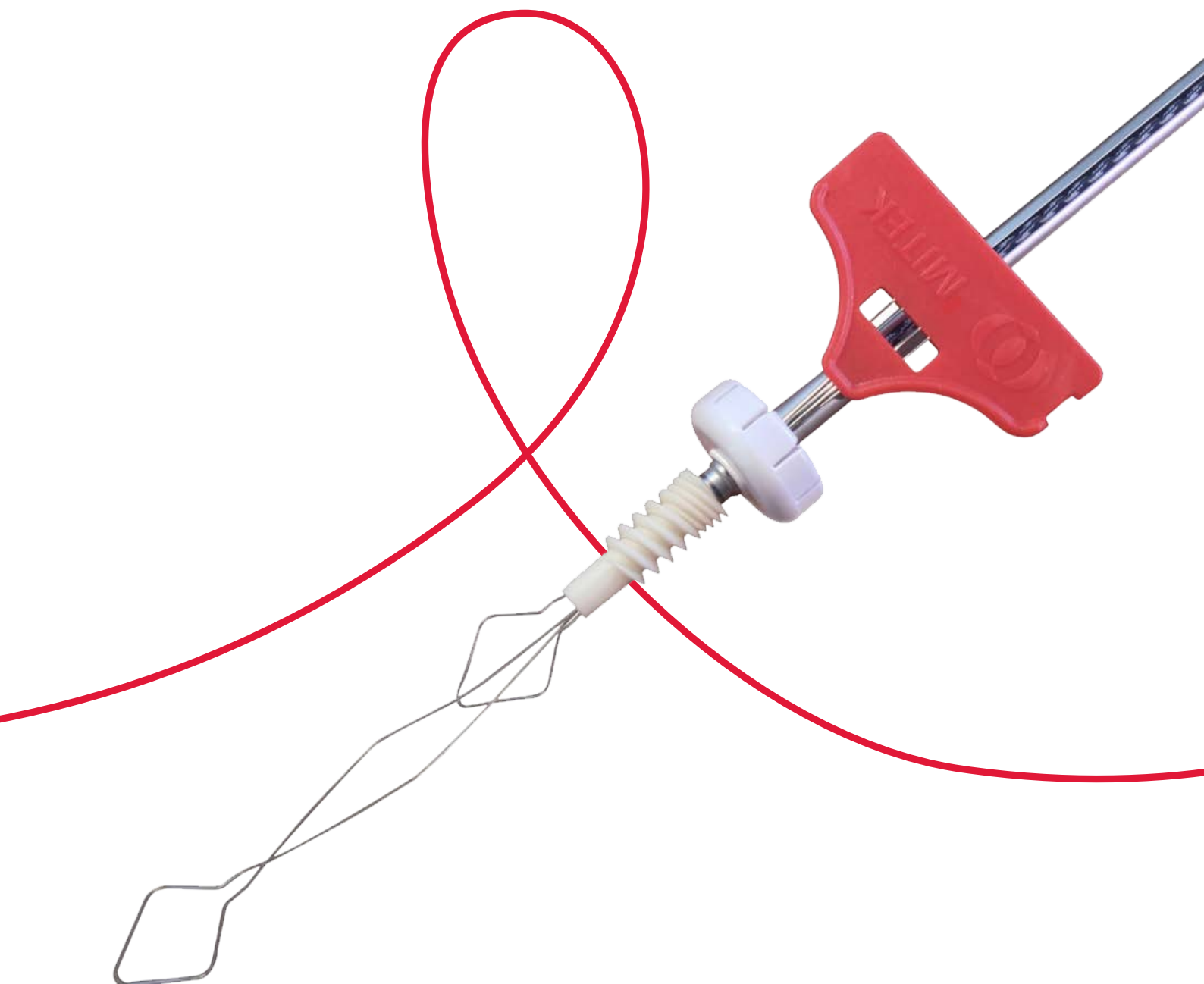

Healix Advance™ Anchors Technology

- Multi-thread Design
- Secure fixation demonstrated through high pullout strength¹
- Cannulated to channel blood flow to the surface`

Versatility

Versatile Suture Configurations

Size & material

- Implants available in 4.75, 5.5, and 6.5 mm sizes.
- Implants available in PEEK and BIOCRYL RAPIDE™ Biocomposite Material.

Dual Wire Threader

Allows increased suture passes through the anchor

- Loads up to 6 strands of #2 ORTHOCORD™ Suture
- 5 strands of #2 DYNACORD™ Suture

Usability

The Sliding Suture Ring

Suture wrapping

- Prevents suture wrapping inside the cannula, while automatically unwinding stitches from above the ring as the inserter is removed.

Suture tensioning

- 360° cleats enable suture management during tensioning.

Stay-Stitch Sliding #2 ORTHOCORD™ Suture

- Can be incorporated into the repair.

Locking suture ring

- For open use



Product Codes

HEALIX ADVANCE Anchors

Item Code	Item Description
222581	4.75mm HEALIX ADVANCE Knotless BR Anchor
222582	5.5mm HEALIX ADVANCE Knotless BR Anchor
222583	6.5mm HEALIX ADVANCE Knotless BR Anchor
222584	4.75mm HEALIX ADVANCE Knotless PEEK Anchor
222585	5.5mm HEALIX ADVANCE Knotless PEEK Anchor
222586	6.5mm HEALIX ADVANCE Knotless PEEK Anchor

HEALIX ADVANCE Instrumentation

Item Code	Item Description
222314	4.5 HEALIX ADVANCE Anchor and 4.75 HEALIX ADVANCE Knotless Anchor Awl
222315	5.5/6.5 HEALIX ADVANCE Anchor Awl
222223	5.5/6.5 HEALIX ADVANCE Knotless Anchor Awl
222316	4.5 HEALIX ADVANCE Anchor Cortical Awl/Tap
222317	5.5/6.5 HEALIX ADVANCE Anchor Cortical Awl/Tap
222328	4.5 HEALIX ADVANCE Anchor Full Tap
222329	5.5/6.5 HEALIX ADVANCE Anchor Full Tap



Purple handle (222314), to be used with 4.75 mm HEALIX ADVANCE Knotless Anchors.



Grey handle (222223), to be used with 5.5 and 6.5 mm HEALIX ADVANCE Knotless Anchors.

References

1. Barber FA, Herbert MA, Beavis RC, Barrera Oro FB. Suture anchor materials, eyelets and designs. update 2008. Arthroscopy. 2008;24(8):859.867.

WARNING: Please also refer to the package insert(s) or other labeling associated with the devices identified in this brochure for additional information.

The third party trademarks used herein are the trademarks of their respective owners.

Please refer to the instructions for use for a complete list of indications, contraindications, warnings and precautions.



Manufactured or distributed by one of more of the following:

Medos International SÀRL
Chemin-Blanc 38
CH-2400 Le Locle
Switzerland

DePuy Mitek, Inc.
325 Paramount Drive
Raynham, MA 02767
T. +1 (800) 382-4682

Johnson & Johnson Medical Pty Ltd
t/a DePuy Synthes
1-5 Khartoum Road,
North Ryde, NSW, 2113
Australia