


Part of an innovative portfolio designed for hernia repair from Ethicon

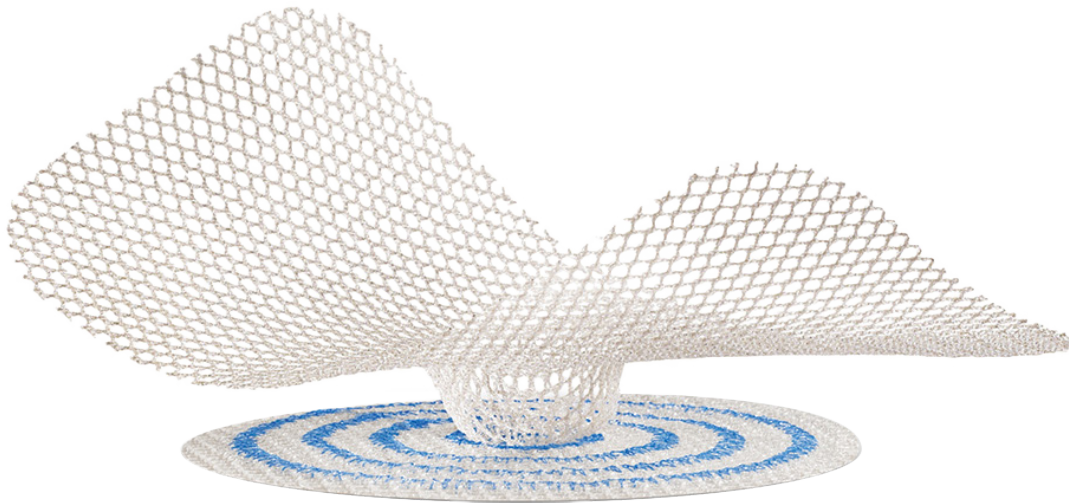

- May offer increased comfort and flexibility¹⁻³
 - Demonstrated significant reduction of pain—both short- and long-term— versus heavyweight mesh^{4*}
- Low rate of recurrence (<1%) at one year postsurgery^{5†}
- Unique MONOCRYL™ (poliglecaprone 25) and PROLENE™ Polypropylene filament construction for enhanced handling and placement^{6,7}
 - MONOCRYL film on underlay facilitates even deployment in the preperitoneal space⁷
- Proven bilayer design combines anterior and posterior repair for long-term security⁸
- ULTRAPRO™ Partially Absorbable Lightweight Mesh leaves behind 65% less foreign material compared to heavyweight mesh after partial absorption⁹

Core mesh material ⁷	ULTRAPRO™ Partially Absorbable Lightweight Mesh (macroporous, thin-filament polypropylene mesh)
Tissue-separating material ⁷	None (Note: This device must always be separated from the abdominal cavity by peritoneum)
Pore size ¹⁰	3.2 mm
Average filament size ¹¹	3.5 mils (pre-absorption of absorbable filament)
Weight ⁹	28 g/m ² (post-absorption)
Customizable ⁷	Yes (patch only)
Fixation requirement ⁷	Adequate onlay patch fixation is required to minimize postoperative complications and recurrence. The fixation technique, method and products used should follow the current standard of care. When fixating with sutures or other mechanical devices, a safe distance from the edge of the mesh of not less than 1 cm must be maintained. The tails created to accommodate the cord structures need to be fixated.
Shelf life	2 years

ULTRAPRO™ Hernia System – Product Specifications

Product Code	Description	Onlay Patch Size	Underlay Patch Size	How Supplied
UHSL	Large	2.4" x 4.7" (6cm x 12cm)	4.0" (~10cm)	Sterile, 3 per box
UHSL6	Large	2.4" x 4.7" (6cm x 12cm)	4.0" (~10cm)	Sterile, 6 per box
UHSM	Medium	2.4" x 4.7" (6cm x 12cm)	3.0" (7.5cm)	Sterile, 3 per box
UHSM6	Medium	2.4" x 4.7" (6cm x 12cm)	3.0" (7.5cm)	Sterile, 6 per box
UHSOV	Oval	2.4" x 4.7" (6cm x 12cm)	4.0" x 4.7" (~10cm x 12cm)	Sterile, 3 per box

For complete indications, contraindications, warnings, precautions, and adverse reactions, please reference full package insert.

References

1. Cobb WS, Kercher KW, Heniford BT. The Argument for Lightweight Polypropylene Mesh in Hernia Repair. *Surg Innov.* 2005;12:63-69. 2. Klosterhalfen B, Junge K, Klinge U. The Lightweight and Large Porous Mesh Concept for Hernia Repair. *Expert Rev Med Devices.* 2005;2:103-17. 3. Post S, Weiss B, Willer M, Neufang T, Lorenz D. Randomized Clinical Trial of Lightweight Composite Mesh for Lichtenstein Inguinal Hernia Repair. *Br J Surg.* 2004;91(1):44-48. 4. Tsirlina VB, Colavita PD, Belyansky I, Zemlyak AY, Lincourt AE, Heniford BT. ULTRAPRO Hernia System Versus Conventional Lichtenstein Repair for Inguinal Hernia System: Results from the Multinational Registry. Abstract presented at 5th International Hernia Congress; March 28-31, 2012; New York, NY. 5. Lorenz R, Koch A, Wiese M, Born H. Gilbert Repair – A New Gold Standard? First Results of 2792 Patients of a Multicentric German Quality Control Study. Poster presented at 15th Annual Hernia Repair; March 13-16, 2013; Orlando, FL. 6. Project PE 02-527. Validation Report: ULTRAPRO Preproduction. June 25, 2003. Ethicon, Inc. 7. ULTRAPRO® Hernia System, Instructions for Use. Ethicon, Inc. 8. Gilbert AI, Young J, Graham MF, Divilio LT, Patel B. Combined Anterior and Posterior Inguinal Hernia Repair: Intermediate Recurrence Rates with Three Groups of Surgeons. *Hernia.* 2004;8(3):203-207. 9. Ethicon Hernia Mesh Specifications. Ethicon, Inc. 10. Cooper S. Test Report: Pore Size ULTRAPRO Annealed. August 3, 2016. Ethicon, Inc. 11. Material Specification for Monofilament Polypropylene, Pigmented or Unpigmented, Regular and Reduced Tip. MS729-002 Rev. 87. June 30, 2015. Ethicon, Inc.

ETHICON

Johnson & Johnson SURGICAL TECHNOLOGIES

©Ethicon US, LLC. 2023. All Rights Reserved. 049286-230622